

BagMixer® 400

400 mL lab blenders

HUBERLAB.

committed to science



interscience

Why choose BagMixer® 400!

Best-seller

- Sold in more than 85 countries
- 35 years experience

High-quality

- Optimal analysis with the sterile **BagFilter®** and **BagPage®**
- In compliance with ISO 7218 standard
- Made in France by **interscience**

Brushless motor

- No wearing parts
- - 50% power consumption
- Extra silent

Safety

- The handle closes the blender and starts the blending
- **Q-Tight®**: no leaks from the top of the bag
- Security drip tray

PATENTED

Click & Clean® removable paddles

- Remove the paddles and easily clean the blending chamber



Easy cleaning

- All stainless steel
- Door opens at 270°
- Removable **Click & Clean®** paddles

Guarantee

- Shock absorbers and window door **life-time** guarantee
- 3 year guarantee (after registration)

PATENTED


Adjustable blending power

- Up to 28 kg of pressure: the hardest samples don't resist! e.g.: nuts, pills, seeds, dog biscuits etc.



Optimal extraction

- From as little as 30 seconds
- Easy control with the digital screen (speed, time, power)

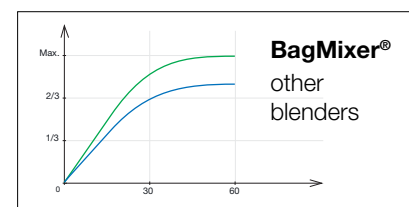

More than 14,000 **BagMixer®** in use every day!

Made in France by **interscience**,
the only manufacturer of blenders and sterile bags for microbiological analysis

Some features are limited to specific models.
For complete technical info, please check online:
www.interscience.fr



BagMixer® 400 models:

The **BagMixer®** lab blender guarantees optimal extraction by the action of its two paddles. Adapted to any sample (food, agricultural, chemical, medical, environmental) to be analyzed between 50 mL and 400 mL, it guarantees homogenization in 30 to 60 seconds.



BagMixer® 400 P
➤ The essential

Ref. 021 230



BagMixer® 400 W
➤ Window door

Ref. 022 230



BagMixer® 400 CC®
➤ Click & Clean®

Ref. 024 230



BagMixer® 400 S
➤ Silent blending

Ref. 025 000



BagMixer® 400 SW
➤ Window door
& LiquidSensor

Ref. 025 100

Also available: packs including bags and accessories; 100 mL and 3500 mL lab blenders. For more info, please contact us.
For more info on BagMixer® S & SW models, please refer to BagMixer® S & SW brochure.

BagMixer® packs:



Starter Pack

- 1 **BagMixer®**
- + **Sterile bags:** 500 BagFilter®
+ 1,000 BagLight®
- + **Accessories:** 50 BagClip®
+ 1 BagOpen® + 1 BagRack®



Classic Pack


- 1 **BagMixer®**
- + **Sterile bags:** 500 BagFilter®
+ 1,000 BagLight® + 500 BagPage®
- + **Accessories:** 50 BagClip®
+ 1 BagOpen® + 1 BagRack®
+ 1 PetriPile® + 1 BagPipet®
+ 1,000 BagTips® Regular
+ 6 Clinet® 101




Premium Pack

- 1 **BagMixer®** & 1 **DiluFlow®**
- + **Sterile bags:** 1,000 BagFilter®
+ 2,000 BagLight® + 1,000 BagPage®
- + **Accessories:** 50 BagClip®
+ 1 BagOpen® + 2 BagRack®
+ 2 PetriPile® + 1 BagPipet®
+ 1,000 BagTips® Regular
+ 6 Clinet® 101

BagSystem® complete range for sample preparation before analysis

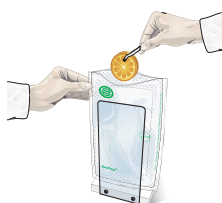

BagFilter®
Collect

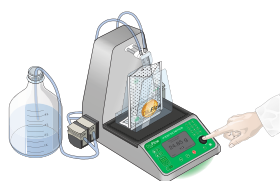

DiluFlow®
Dilute








BagMixer®
Blend




BagTools®
Pipet

Technical specifications

		BagMixer® 400 P	BagMixer® 400 W	BagMixer® 400 CC®
	Reference	021 230	022 230	024 230
Key features	Removable and autoclavable paddles, Click & Clean® system	-	-	✓
	Blending power indicator	-	-	✓
	Adjustable blending power (up to 28 kg of pressure)	-	-	✓
	Digital screen	-	-	✓
	Leak free Q-tight® system	-	-	✓
	Side-by-side paddle stop	-	-	✓
	Silent brushless motor	-	-	✓
	Window door	-	✓	✓
	Security drip tray	in option	✓	✓
	270° opening door	✓	✓	✓
	All stainless steel	✓	✓	✓
	Compact and ergonomic	✓	✓	✓
Technical specifications	Door weight	1050 g	950 g	
	Blending time	30, 60, 90, 120, 150, 180, 210 sec., ∞		from 1 sec. to 59 min, ∞
	Blending speed	fixed speed: 8 strokes/sec.		variable speed: 4, 6, 8, 10 strokes/sec.
	Dimension (w x d x h)	26 x 39 x 29 cm		26.5 x 42 x 26 cm
	Weight	16.5 kg		16 kg
	Power	~110/220 ~ 50-60Hz		
	Capacity	50-400 mL		
Guarantees	Made in France	✓	✓	✓
	Life-time guarantee	Unbreakable shock absorbers and window door		
	3 year guarantee (after registration)	✓	✓	✓
Manufactured in compliance with				  

Discover our complete range for microbiology

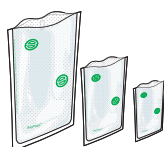

BagFilter®
Sterile bags




DiluFlow®
Gravimetric dilutors



FlexiPump®
Peristaltic pumps



easySpiral Dilute®
Automatic diluter et platers

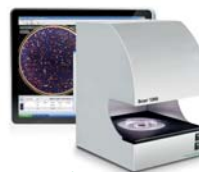

Scan®
Automatic colony counters