

**DURAN®**


STAINLESS STEEL

1-, 2- OR 3-PORT

CONNECTOR CAP

GL 45

committed to  
**brands.**

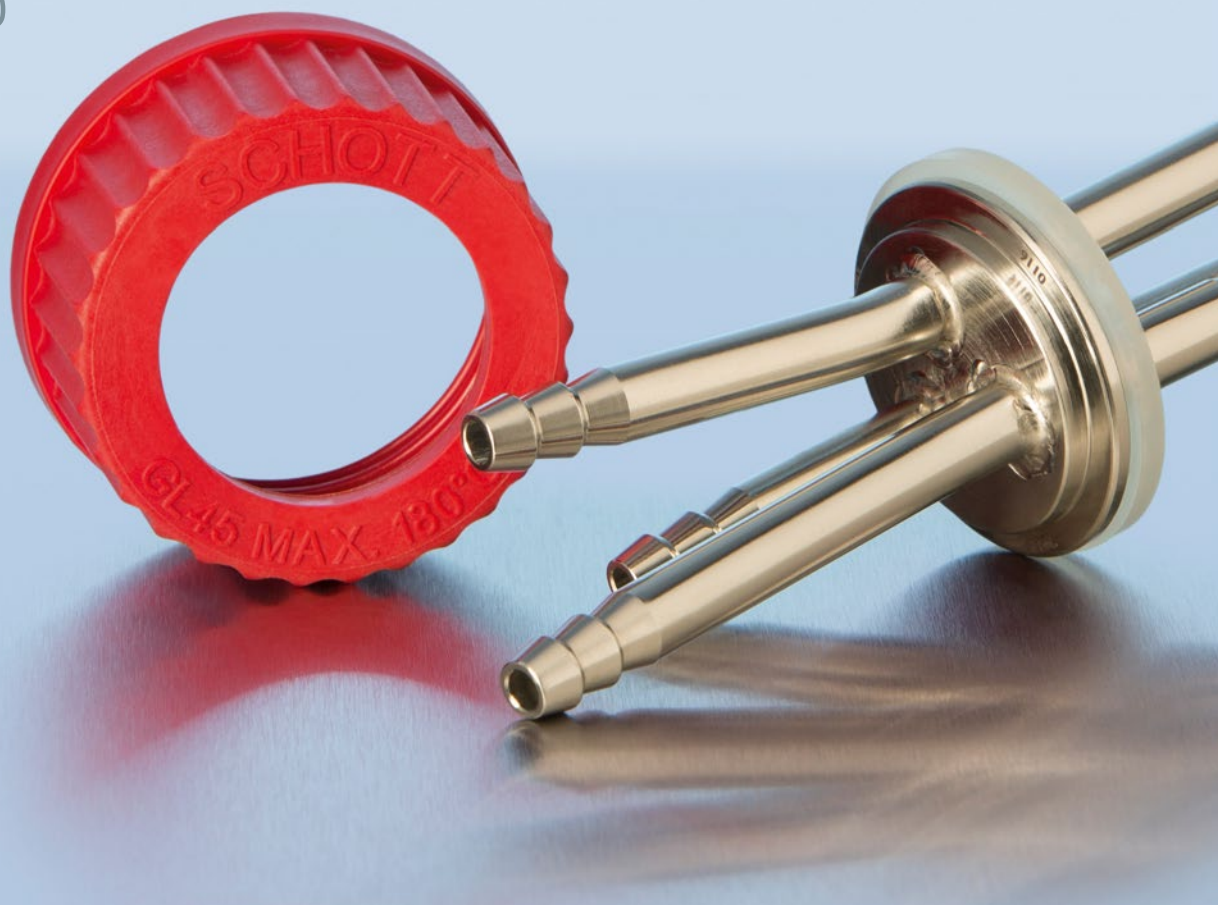

**HUBERLAB.**

committed to science



**DURAN  
WHEATON  
KIMBLE**

Excellence in your hands

# DURAN® STAINLESS STEEL 1-, 2- OR 3-PORT CONNECTOR CAP GL 45

A robust and durable tubing connection system for use with all sizes of DURAN® GL 45 laboratory and media bottles. The bottle top adaptor facilitates the transfer of liquids within a closed and sterile system. Ideal for chemistry, life science and biopharma laboratories.

## Applications

- Sterile liquid transfer between containers
- Liquid sampling into bottles from larger containers
- Reservoir cap for supplying laboratory instruments
- Port assemblies for spinner flasks, and cell culture vessels
- Autoclaving media in GL 45 DURAN® bottles



## Features

- GL 45 connector cap can be used for all sizes of DURAN® GL 45 bottles from 100 to 20,000 ml
- The one, two or three barbed, male hose connectors on both the top and bottom assemblies are suitable for flexible tubing with a 8.0 mm inner diameter, such as Sani-Tech® Ultra (ULTRA-C-312) Silicone tubing\*
- The 1-port cap is suitable for autoclaving media when used with a sterile air venting capsule filter
- Freely rotating stainless steel centre allows the bottle to be unscrewed without the need to disconnect the tubing
- The two or three top connectors are slightly splayed to facilitate easier attachment of the tubing
- Connectors and body are manufactured from Type 316L (1.4404 / S31603) stainless steel
- All components are fully welded to avoid entrapment voids
- Hygienic, electropolished surface for efficient cleaning and sterilisation
- A liquid tight seal is provided by a (Platinum catalyzed USP <88> Class VI) silicone gasket and the red PBT GL 45 screw cap
- Wide temperature usage range (up to 180 °C)
- Autoclavable (@ 121 °C / 15 minutes)
- Retrace coded for full traceability

\*Sani-Tech® is a registered trademark of Saint-Gobain

## Ordering information

Description	Pack Unit	Cat. No.
<b>NEW</b> DURAN® Stainless steel 1-port connector cap GL 45	1	29 261 26 07
DURAN® Stainless steel 2-port connector cap GL 45	1	29 261 27 01
DURAN® Stainless steel 3-port connector cap GL 45	1	29 261 28 04

## Replacement parts

Description	Pack Unit	Cat. No.
DURAN® Silicone Sealing Gasket 40.5 mm outer and 27.5 mm inner diameter, 3 mm thick	10	29 223 28 05
DURAN® Screw Cap GL 45, red, PBT with 34 mm central aperture	10	29 227 10 07

These products are available from your specialist laboratory dealer.



# Contact.

## HUBERLAB.

committed to science

**OneStop. OneShop – HUBERLAB. bietet das Rundum-Sorglos-Paket für alle Bedürfnisse im Labor.**

Profitieren Sie von unseren langjährigen Partnerschaften zu Top-Lieferanten.

Als Vollversorger für das Labor lässt unser umfassendes Sortiment keine Wünsche offen.

---

**OneStop. OneShop – HUBERLAB. Nous offre le package sans souci pour tous les besoins de votre laboratoire.**

Profitez de nos partenariats de longue date avec les meilleurs fournisseurs.

En tant que fournisseur complet pour le laboratoire, notre gamme complète ne laisse rien à désirer.

---

**OneStop. OneShop – HUBERLAB. offers the all-included service package for all of your laboratory requirements.**

Benefit from our long-standing partnerships with top-suppliers.

As a full-range supplier for the laboratory, our comprehensive range leaves nothing to be desired.



**Wir beraten Sie gerne persönlich**  
Nous vous conseillons personnellement  
We'll be happy to advise you personally



T +41 61 717 99 77



F +41 61 711 93 42



info.huberlab.ch



www@huberlab.ch



Let's chat  
www.huberlab.ch



**HUBERLAB. AG**  
Industriestrasse 123  
4147 Aesch