

committed to
brands.



**DURAN
WHEATON
KIMBLE**

Excellence in your hands

PART OF THE **DURAN®** BOTTLE SYSTEM

NEW: DURAN® GL 45 CONNECTION CAP WITH 2, 3, 4, OR 5 STAINLESS STEEL CONNECTORS FOR NARROW TUBING

The reusable DURAN® GL 45 Connection Cap is ideal for **aseptic fluid transfer** into or out of DURAN® bottles using peristaltic pump tubing. The screw caps will fit all DURAN® Original GL 45 laboratory bottles 100 mL and larger. The barbed tubing connectors (4.2 mm outer diameter) accept **narrow BPT peristaltic pump tubing** with a 3.2 mm (1/8") inner diameter.

FEATURES & BENEFITS

- **Versatile** choice of 2, 3, 4, or 5 stainless steel (316L) **tubing connectors**
- **Flexibility** with connectors for internal drop tubes
- Manufactured from **hygienic, chemically resistant materials**
- **Sterilizable / autoclavable** at +121 °C for 20 minutes
- **High temperature resistance** up to +180 °C
- **Reliable leak-resistant sealing** with fluoroelastomer O-Rings
- Construction materials certified **as food contact safe**
- **Tangle-free tubes** with the PTFE center freely rotates
- Robust and reusable with a **long service life**



PERISTALTIC PUMPING – WHAT IS THAT?


Peristaltic pumps use flexible tubing that run through rollers in the pump head. Peristaltic pumps confine the media inside the tubing, so that the liquid cannot contaminate the pump, and the pump cannot contaminate the liquid. This means that this method of pumping liquids is particularly suitable for use in high purity applications. The gentle, low shear pumping action is ideal for handling suspensions of cells or large biomolecules in the laboratory.

NEW: DURAN® GL 45 CONNECTION CAP

WITH 2, 3, 4, OR 5 STAINLESS STEEL CONNECTORS FOR NARROW TUBING

APPLICATION AREAS

- Sample collection
- Bioreactor feeds
- Buffer and media addition
- Peristaltic pumping of fluids
- Perfusion circuit for cell culture
- Use in a chemistry micro-reactor
- Oxygenation or aeration of samples

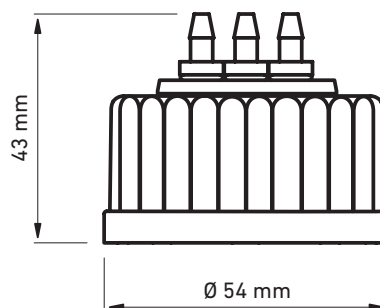

ORDERING INFORMATION

Description	Pack Unit	Cat. No.
DURAN® GL 45 Connection Cap with red PBT screw cap, PTFE insert with 2-connectors (stainless steel)	1	29 263 28 06
DURAN® GL 45 Connection Cap with red PBT screw cap, PTFE insert with 3-connectors (stainless steel)	1	29 264 28 07
DURAN® GL 45 Connection Cap with red PBT screw cap, PTFE insert with 4-connectors (stainless steel)	1	29 265 28 08
DURAN® GL 45 Connection Cap with red PBT screw cap, PTFE insert with 5-connectors (stainless steel)	1	29 266 28 09



Discover our wide range of borosilicate glass 3.3 laboratory bottles in different sizes, with closures, liquid transfer caps and accessories.



Explore more:
www.duran-bottle-system.com



Contact.

HUBERLAB.

committed to science

OneStop. OneShop – HUBERLAB. bietet das Rundum-Sorglos-Paket für alle Bedürfnisse im Labor.

Profitieren Sie von unseren langjährigen Partnerschaften zu Top-Lieferanten.

Als Vollversorger für das Labor lässt unser umfassendes Sortiment keine Wünsche offen.

OneStop. OneShop – HUBERLAB. Nous offre le package sans souci pour tous les besoins de votre laboratoire.

Profitez de nos partenariats de longue date avec les meilleurs fournisseurs.

En tant que fournisseur complet pour le laboratoire, notre gamme complète ne laisse rien à désirer.

OneStop. OneShop – HUBERLAB. offers the all-included service package for all of your laboratory requirements.

Benefit from our long-standing partnerships with top-suppliers.

As a full-range supplier for the laboratory, our comprehensive range leaves nothing to be desired.



Wir beraten Sie gerne persönlich
Nous vous conseillons personnellement
We'll be happy to advise you personally



T +41 61 717 99 77



F +41 61 711 93 42



info.huberlab.ch



www@huberlab.ch



Let's chat
www.huberlab.ch



HUBERLAB. AG
Industriestrasse 123
4147 Aesch